


# EHCA1200 system

## Cellogel cellulose acetate system

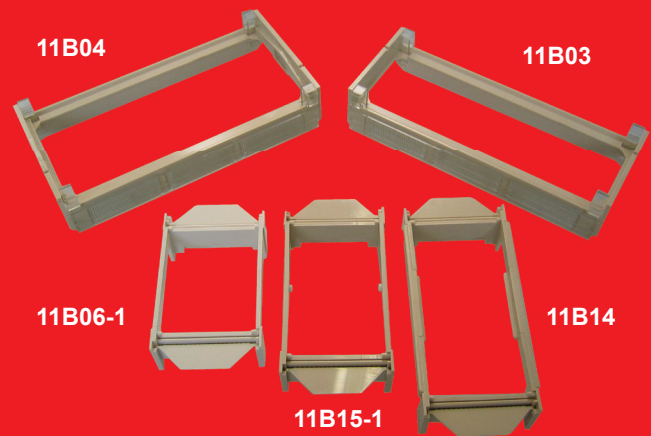

Cellogel is a film of cellulose acetate in gel form. Cellogel is the ideal electrophoretic support for clinical electrophoresis and for the immunological techniques.

- Compact high resolution system for clinical electrophoresis
- Designed for routine and research needs
- Easy loading with bridges
- Fully compatible with Cellogel precast gels and kits
- Complete range of cellulose acetate gels and kits
- Densitometer software and scanner available

### Recommended power supply


EV0220  
EV1450

Functions with six strips 2.5x14 cm or with three strips 5.7x14 cm on three bridges, model France of 8.5 cm, furthermore it works with the same size strips as above and with sheets 14x14 cm or 18.3x14 cm on an 8.5 cm long bridge with plastic clips.

Injection moulded polycarbonate with high chemical and physical resistance. The lid is in semi-transparent polycarbonate with two magnets which work safety micro-switches and cut off the current when the lid is taken off.

Universal bridges supports Cellogel during sample application by serving as a convenient loading template for the required applicator. Bridges also available for Helena-, Pratiga- and Shandon-type strips and other size formats.

A wide selection of precast gel and strips are available.

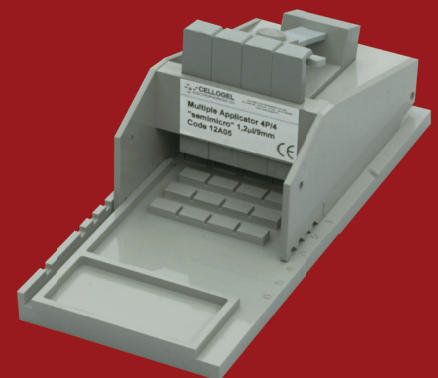

Code	Description
EHCA1200-SYS	Horizontal unit for cellulose acetate including 3 bridges 8.5 for Cellogel strips 2.5 x 14 cm and 5.7 x 14 cm, 1 bridge for Cellogel sheets 18.3 x 14 cm
EHCA1200-BR11B06-1	France bridge, 8.5 cm, for strips of 2.5x14 cm or 5.7x14 cm
EHCA1200-BR11B03	Long bridge for sheets of 18.3x14 cm or 14x14 cm and strips of 2.5x14 cm or 5.7x14 cm. Migration field 8.5 cm
EHCA1200-BR11B15-1	France bridge, 11 cm, for strips of 2.5x17 cm or 5.7x17 cm, Rectangular.
EHCA1200-BR11B04	Long bridge for sheets of 17x17 cm and strips of 5.7x17 cm. Migration field 11 cm.
EHCA1200-BR11B14	France bridge, 14 cm, for Cellogel RS Wedge of 5x18.5 cm and 5.7x18.5 cm Rectangular.
EHCA1200-AP08-8P2	8 Sample Micro Applicator
EHCA1200-AP02-SU	2+2 Samples Semimicro Applicator for 2/5 x 14 cm and IFE kit
EHCA1200-AP08-4P4	4 Samples Semimicro Applicator
EHCA1200-AP05	4 Samples Semimicro Applicator
EHCA1200-AP08-4CS	4 Samples Semimicro Applicator for USP CHONDROITIN SULFATE test
EHCA1200-AP08-6P2	6 Samples Semimicro Applicator
EHCA1200-KC30-R	Serum Proteins kit
EHCA1200-KC31	High Resolution Serum Proteins kit
EHCA1200-KC09	IEF Serum + Concentrated urine kit
EHCA1200-KC64	Glycosylated Hemoglobins HbA1c kit
EHCA1200-KC09-02	Immunofixation
EHCA1200-KC35	Hemoglobins
EHCA1200-KC42	Lipoproteins
EHCA1200-SOFT	Turboscan Universal Densitometer Software
EHCA1200-SCAN	Scanner for EHCA1200-SOFT



Artistic in science

Consort bvba  
Parklaan 36  
B-2300 Turnhout  
Belgium

Telephone [\(+32\) \(0\)14 41 12 79](tel:+322014411279)  
Email [info@consort.be](mailto:info@consort.be)  
Sales [sales@consort.be](mailto:sales@consort.be)  
Support [support@consort.be](mailto:support@consort.be)  
Website [www.consort.be](http://www.consort.be)