


- Quick & Efficient
- Convenient & Versatile
- Safe & Reliable

Quick transfer times make these systems ideal for the rapid and efficient transfer of nucleic acids and proteins from agarose or acrylamide gels.

Solid plate style electrodes assure even pressure and complete molecular transfer.

Galileo semi-dry electroblotters will accommodate most pre-cast gels, and for many larger format pre-cast gels, once the stacking area of the gel is removed the "working" area will also fit easily into the devices. Since only the transfer "sandwich" of gel, membrane and blotting papers must be kept wet, much less buffer is required than traditional tank transfer systems.

Interlocking safety lid prevents reverse orientation of electric field. High quality plate electrodes (stainless steel cathode and platinum anode) and gold plated components assure long, trouble-free performance. System is supplied complete with power leads and a sample pack of our ultra-pure cotton fiber filter paper.

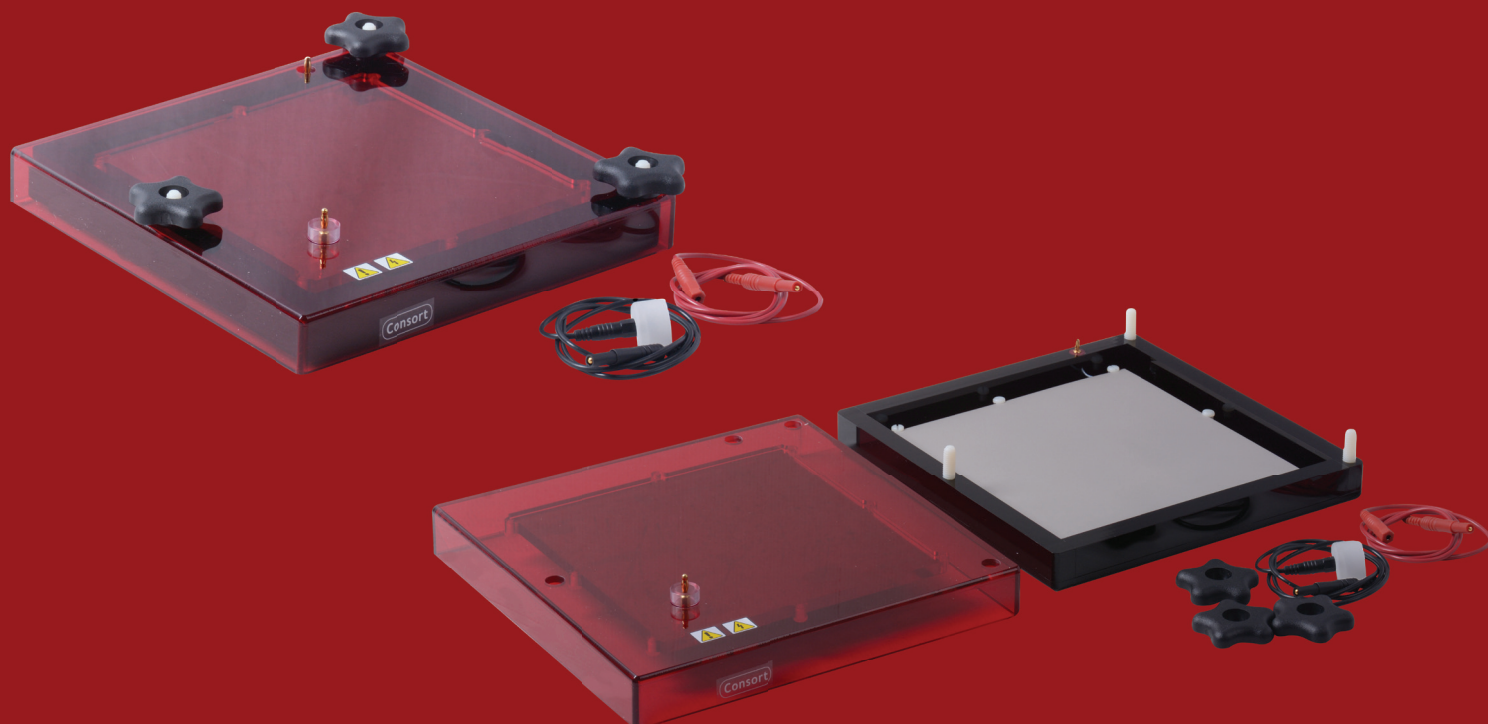

Code	Description
ESDB3100	Semi Dry Electroblotter, 11x11cm blotting area
ESDB3300	Semi Dry Electroblotter, 21x21cm blotting area
EV2310	Power supply, 300 V, 1000 mA, 150 W
EV3020	Power supply, 300 V, 2000 mA, 300 W



Artistic in science

Consort bvba  
Parklaan 36  
B-2300 Turnhout  
Belgium

Telephone [\[+32\] \[0\]14 41 12 79](tel:+322014411279)  
Email [info@consort.be](mailto:info@consort.be)  
Sales [sales@consort.be](mailto:sales@consort.be)  
Support [support@consort.be](mailto:support@consort.be)  
Website [www.consort.be](http://www.consort.be)